



ChibaSono & Chiba

*Fine Needle
Aspiration Biopsy*

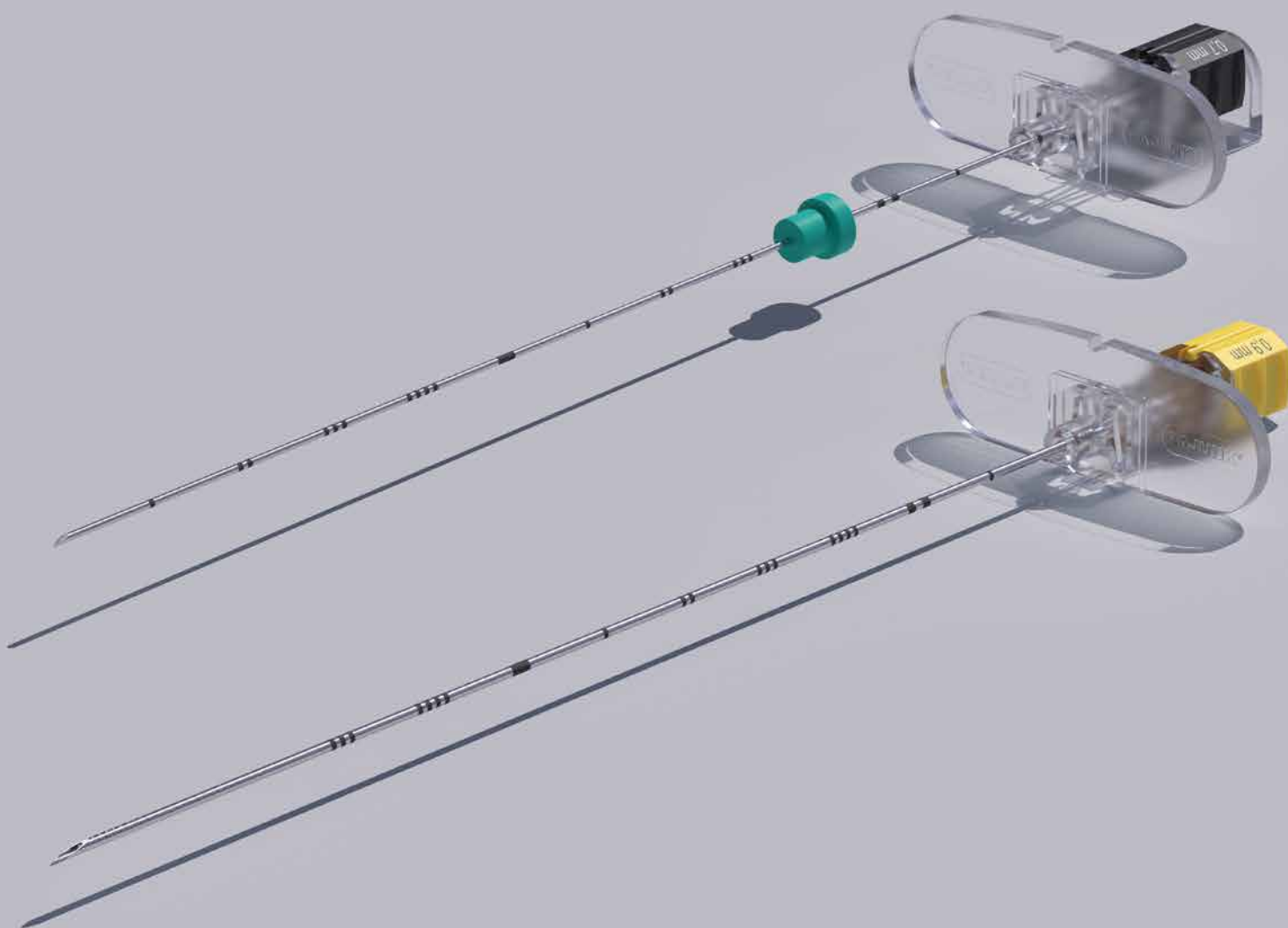

Echogenic & Non-echogenic Chiba Needles

Ultrasound guided biopsies have been increasing in the recent years. Together with Dr. Chris Mitchell, PAJUNK® has developed the ChibaSono, a new generation of needles that allow for a precise localisation under ultrasound. Needle and tip design allow the user to focus on the procedure and patient by making the tip visualisation under ultrasound faster, easier and safer. The wide range of ChibaSono and Chiba needles, sizes and lengths offers the optimal tool for percutaneous cholangiography, cytological biopsy and general percutaneous applications.

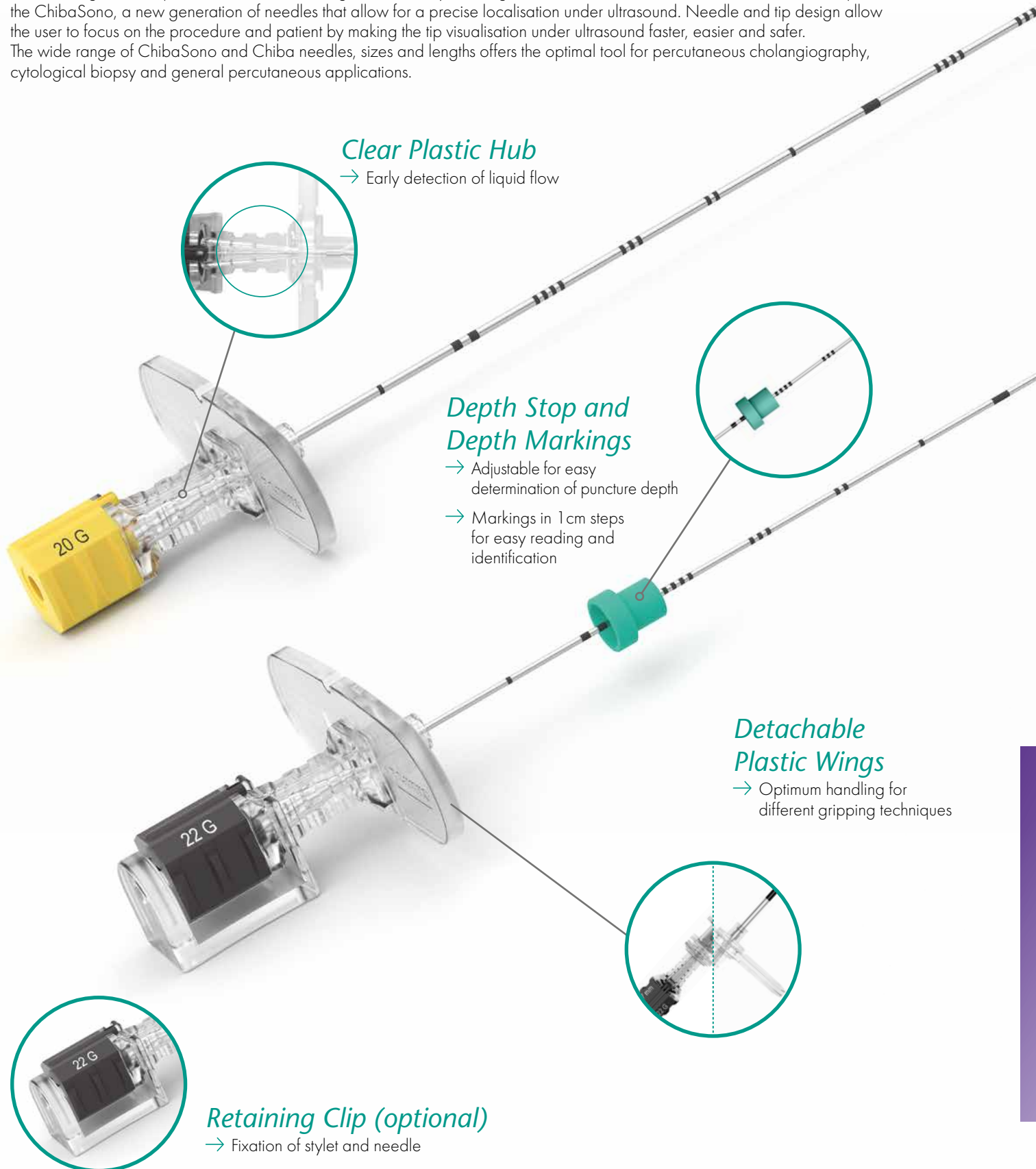

Clear Plastic Hub

→ Early detection of liquid flow

Depth Stop and Depth Markings


- Adjustable for easy determination of puncture depth
- Markings in 1 cm steps for easy reading and identification

Detachable Plastic Wings

→ Optimum handling for different gripping techniques

Retaining Clip (optional)

→ Fixation of stylet and needle

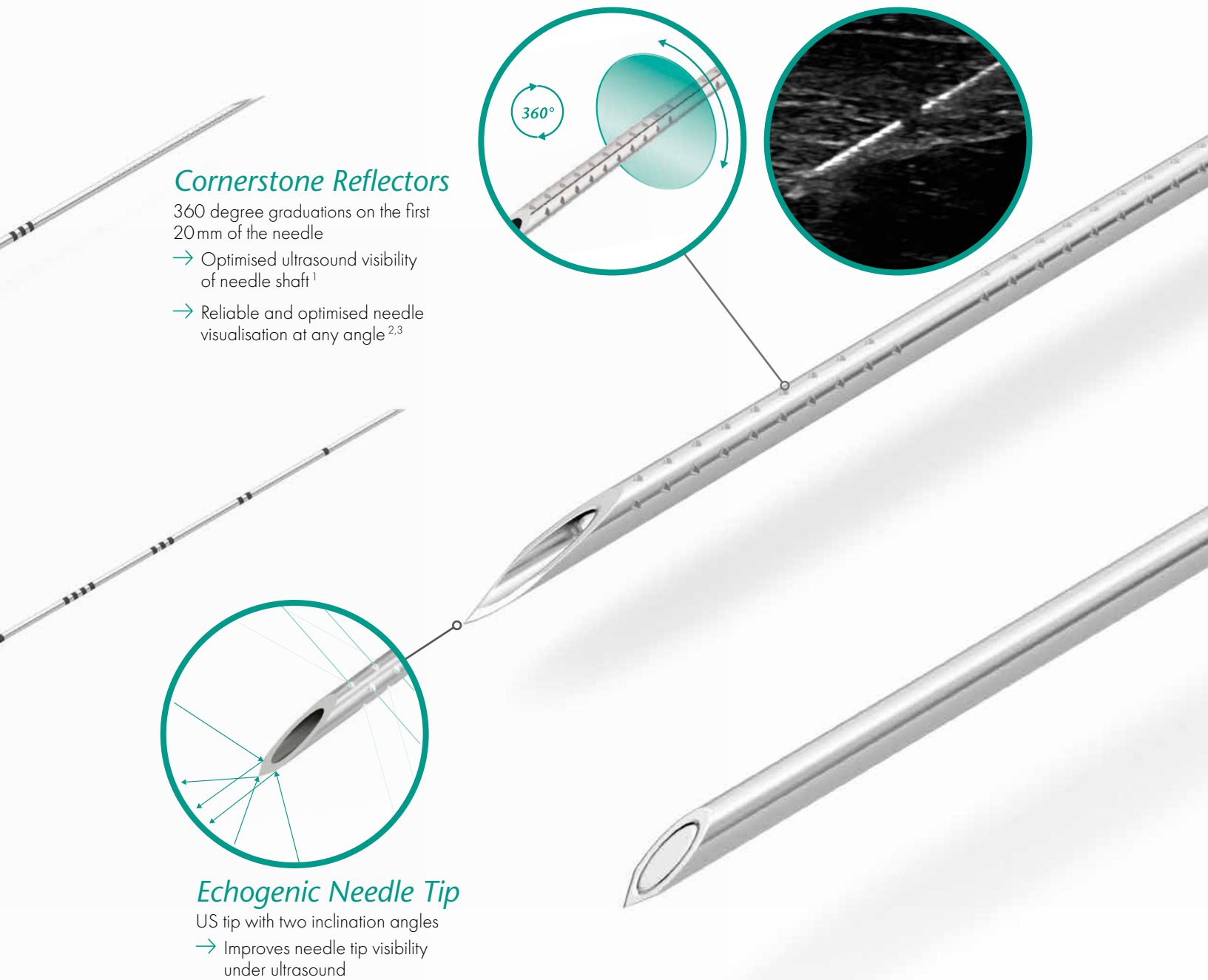

Cornerstone Reflectors

360 degree graduations on the first 20mm of the needle

- Optimised ultrasound visibility of needle shaft¹
- Reliable and optimised needle visualisation at any angle^{2,3}

Echogenic Needle Tip

US tip with two inclination angles

- Improves needle tip visibility under ultrasound



AVAILABLE OPTIONS

Tip options for Chiba:

Beveled Tip

- ▶ Designed for fluoroscopy guided procedures
- ▶ Precision-ground metal stylet
- ▶ Minimized tissue carry-out into the puncture
- ▶ Increased needle stability



US Tip

- ▶ Special grinding with two inclination angles improves needle tip visibility under ultrasound
- ▶ Rounded, blunt stylet – set back from tip



ChibaSono only with US tip

ChibaSono & Chiba

ChibaSono

with US tip, Cornerstone Reflectors, depth markings and wings

Size	Item no.	PU
22G x 50 mm (2")	530S070050	25
22G x 90 mm (3 1/2")	530S070090	25
22G x 100 mm (4")	530S070100	25
22G x 125 mm (5")	530S070125	10
22G x 150 mm (6")	530S070150	10
22G x 220 mm (8 2/3")	530S070220	10
20G x 90 mm (3 1/2")	530S095090	25
20G x 120 mm (4 3/4")	530S095120	10
20G x 150 mm (6")	530S095150	10
20G x 220 mm (8 2/3")	530S095220	10
18G x 150 mm (6")	530S120150	10

Chiba

with beveled tip, retaining clip and wings

Size	Item no.	PU
22G x 90 mm (3 1/2")	201S070090	25
22G x 120 mm (4 3/4")	201S070120	10
22G x 150 mm (6")	201S070150	10

Chiba

with beveled tip, depth stop, depth markings, retaining clip and wings

Size	Item no.	PU
22G x 120 mm (4 3/4")	202S070120	10
22G x 150 mm (6")	202S070150	10
22G x 220 mm (8 2/3")	202S070220	10
20G x 120 mm (4 3/4")	202S095120	10
20G x 150 mm (6")	202S095150	10
20G x 220 mm (8 2/3")	202S095220	10
20G x 280 mm (11")	202S095280	10

Chiba

with US tip, retaining clip and wings

Size	Item no.	PU
22G x 220 mm (8 2/3")	200S070220	10



Scan now!

PAJUNK® GmbH Medizintechnologie
Global Headquarters
Karl-Hall-Straße 1
78187 Geisingen, Germany
Tel.: +49 7704 9291-0
Fax: +49 7704 9291-600
info@pajunk.com
pajunk.com

PAJUNK® Medical Produkte GmbH
Central Europe
Karl-Hall-Straße 1
78187 Geisingen, Germany
Tel.: +49 7704 8008-0
Fax: +49 7704 8008-150
info@pajunk-medipro.com
pajunk.de & pajunk.eu

PAJUNK® Swiss AG
Switzerland
Neuhofstrasse 5A
6340 Baar, Switzerland
Tel.: +41 417 609 627
info@pajunk.ch
pajunk.ch

PAJUNK® Medical System L.P.
United States of America
4575 Marconi Drive
Alpharetta, GA 30005, USA
Tel.: +1 888 972 5865
Fax: +1 678 514 3388
info@pajunk-usa.com
pajunkusa.com

PAJUNK® UK Medical Products Ltd
United Kingdom
Unit D1 The Waterfront
Goldcrest Way, Newburn Riverside
Newcastle upon Tyne, NE15 8NY, UK
Tel.: +44 191 264 7333
info@pajunk.co.uk
pajunk.co.uk