

PAJUNK®

*Dental Syringes
Systems for
dental anaesthesia*



Dental

MADE IN GERMANY

Cylinder cartridge syringes

Dental syringes by PAJUNK®

PAJUNK® offers the dental cylinder cartridge syringes Quickject in different variants: as compact slip-on device or folding mechanism with rotatable finger rest or as a compact version. They are all equipped with a swiveling cylinder cartridge holder. The stainless-steel version Varioject is also available in two compact versions: with slip-on device or folding mechanism. Varioject is manufactured completely of stainless steel. The robust material is particularly resistive against damaging and harsh and aggressive cleaning agents. In addition, the special matted surface finish avoids irritating reflections. The ergonomic design of the syringes permits a safe and convenient handling.

- ➔ Ergonomic design – safe and convenient handling
- ➔ Slide-on technology permits fast replacement during treatment

Varioject made of stainless steel

- ➔ Resistive against damaging and harsh and aggressive cleaning agents
- ➔ Matted surface finish avoids irritating reflections

Whether slip-on device or folding mechanism, both techniques support either manual aspiration by means of three extensible aspiration hooklets or arrowhead as well as self-aspiration. They are suitable for the use of regular commercially available 1.8ml cylinder cartridges, and for cannulas equipped with a metric thread.

Syringes with slip-on device

If using the slip-on technology, several filled cartridge holders can be prepared and fitted.

➔ This guarantees fast exchange during treatment.



Syringes with folding mechanism

In contrast to the slip-on device, the folding technology permits the cartridge holder to be swung sideways, and only the cartridge itself is replaced. Alternatively to the folding mechanism with rotatable finger rest, a compact variant with rigid finger rest.



Folding technology – insertion of the cartridge



Dental syringes

All at a glance

Quickject


Aspiration hooklets



Arrowhead



Cylindrical pressure plate

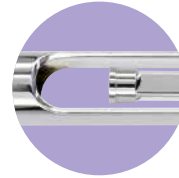

- for cylinder cartridges 1.8 ml
- with thumb ring
- rotatable finger rest

Folding technology, rotatable

| Product | Item no. | PU | Thumb ring |
|---|------------|----|------------|
| With 3 aspiration hooklets | 8990-10004 | 1 | |
| With 3 aspiration hooklets and self-aspiration | 8990-10006 | 1 | |
| With arrowhead | 8990-10005 | 1 | |
| With arrowhead and self-aspiration | 8990-10007 | 1 | |
| With cylindrical pressure plate and self-aspiration | 8990-10002 | 1 | |



- for cylinder cartridges 1.8 ml
- with thumb ring
- with rigid finger rest

Folding technology, compact

| Product | Item no. | PU | Thumb ring |
|---|------------|----|------------|
| With 3 aspiration hooklets | 8990-25004 | 1 | |
| With 3 aspiration hooklets and self-aspiration | 8990-25006 | 1 | |
| With cylindrical pressure plate and self-aspiration | 8990-25002 | 1 | |



- for cylinder cartridges 1.8 ml
- with thumb ring
- with rigid finger rest

Slip-on device, compact

| Product | Item no. | PU | Thumb ring |
|--|------------|----|------------|
| With 3 aspiration hooklets | 8990-20004 | 1 | |
| With 3 aspiration hooklets and self-aspiration | 8990-20006 | 1 | |
| With arrowhead | 8990-20005 | 1 | |
| With arrowhead and self-aspiration | 8990-20007 | 1 | |
| Slip-on cartridge holder for cylinder cartridges 1.8 ml | 8990-00002 | 3 | |
| Slip-on cartridge holder for cylinder cartridges 1.8 ml with self-aspiration | 8990-00012 | 3 | |



Slip-on cartridge holder

Dental syringes

All at a glance

Varioject

Aspiration hooklets



Arrowhead



- made of stainless steel
- for cylinder cartridges 1.8 ml
- with thumb ring
- rigid finger rest



- made of stainless steel
- for cylinder cartridges 1.8 ml
- with thumb ring
- rigid finger rest

Folding technology, compact

Thumb ring

| Product | Item no. | PU |
|--|------------|----|
| With 3 aspiration hooklets and self-aspiration | 8990-50006 | 1 |
| With arrowhead and self-aspiration | 8990-50007 | 1 |

Slip-on device, compact

Thumb ring

| Product | Item no. | PU |
|--|------------|----|
| With 3 aspiration hooklets and self-aspiration | 8990-60006 | 1 |
| Slip-on cartridge holder for Varioject 1.8 ml with self-aspiration | 8990-01012 | 3 |

Components



With round groove, chromium-plated, slitted

Chromium-plated, non-slitted

| Product | Item no. | PU |
|---|------------|----|
| Ointment cannulas, sterile | | |
| Ø 1.0 mm, with detent notch and cleaning stylet | 8990-05000 | 10 |
| Attachment unit | | |
| With round groove, chromium-plated, slitted | 8990-02000 | 1 |
| Chromium-plated, non-slitted | 8990-02010 | 1 |

PAJUNK® GmbH
 Medizintechnologie
 Karl-Hall-Strasse 1
 D-78187 Geisingen / Germany
 Phone +49 (0) 77 04/92 91-0
 Telefax +49 (0) 77 04/92 91-6 00
www.pajunk.com